Development of an Aerosol Concentration Model to Vary Exposure Duration with Infectious Agents

Justin M. Hartings, Ph.D.

Department of Aerobiology & Product Evaluation
United States Medical Research Institute for Infectious Diseases
(USAMRIID)

Fort Detrick, Maryland U.S.A.

19 NOVEMBER 2003

maintaining the data needed, and c including suggestions for reducing	lection of information is estimated to completing and reviewing the collect this burden, to Washington Headqu uld be aware that notwithstanding an DMB control number.	ion of information. Send comments arters Services, Directorate for Infor	regarding this burden estimate mation Operations and Reports	or any other aspect of the s, 1215 Jefferson Davis	nis collection of information, Highway, Suite 1204, Arlington	
		2. REPORT TYPE N/A		3. DATES COVERED		
4. TITLE AND SUBTITLE				5a. CONTRACT NUMBER		
Development of an Aerosol Concentration Model to Vary Exposure Duration with Infectious Agents				5b. GRANT NUMBER		
Duration with infectious Agents				5c. PROGRAM ELEMENT NUMBER		
6. AUTHOR(S)				5d. PROJECT NUMBER		
				5e. TASK NUMBER		
				5f. WORK UNIT NUMBER		
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Department of Aerobiology & Product Evaluation United States Medical Research Institute for Infectious Diseases (USAMRIID) Fort Detrick, Maryland U.S.A. 8. PERFORMING ORGANIZATION REPORT NUMBER						
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)		
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)		
12. DISTRIBUTION/AVAIL Approved for publ	LABILITY STATEMENT ic release, distributi	on unlimited				
	otes 51, Proceedings of t Research, 17-20 No					
14. ABSTRACT						
15. SUBJECT TERMS						
16. SECURITY CLASSIFICATION OF: 17. LIMITATION				18. NUMBER	19a. NAME OF	
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified	ABSTRACT UU	OF PAGES 31	RESPONSIBLE PERSON	

Report Documentation Page


Form Approved OMB No. 0704-0188

Aerobiology Mission

- Develop appropriate animal models
 - Define pathogenesis/mechanism of toxicity/toxicology
 - Develop surrogate markers of efficacy
- Develop new bioaerosol technologies
 - Precision and accuracy of dose control
 - Aerosol size
 - Environmental conditions

Aerosol Challenge

- Deliver challenge agent by aerosol
- Under desired conditions
 - Temperature
 - Humidity
 - Particle size
- At the required dose

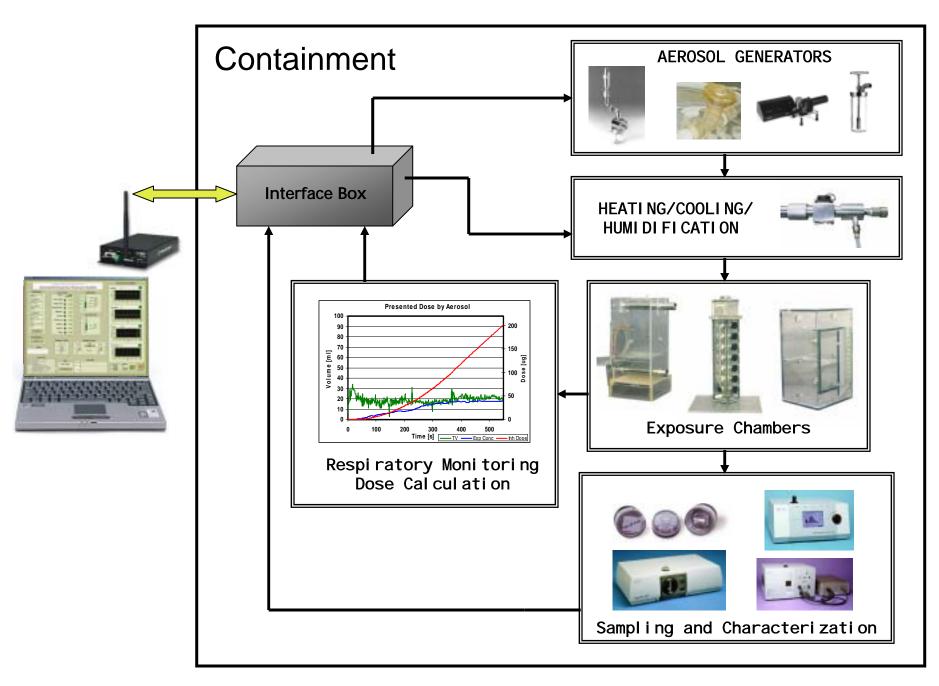

Integrated Aerosol Control and Management Platform

System Hardware



Workstation



Containment

Challenge Dose

Aerosol Concentration

Respiratory
Function

Exposure Duration

Challenge Dose

$$D(t_{\exp}) = \int_0^{t_{\exp}} R(t)C(t)dt$$


Toxin Dose Control

• Calculate dose integral in real-time

$$D(t_{\exp}) = \int_0^{t_{\exp}} R(t)C(t)dt$$

- Exposure duration dictated by measurements
- Automatic compensation for variations in R(t) and C(t)
- Real-time dose calculation controls t_{exp}

Chamber Concentration

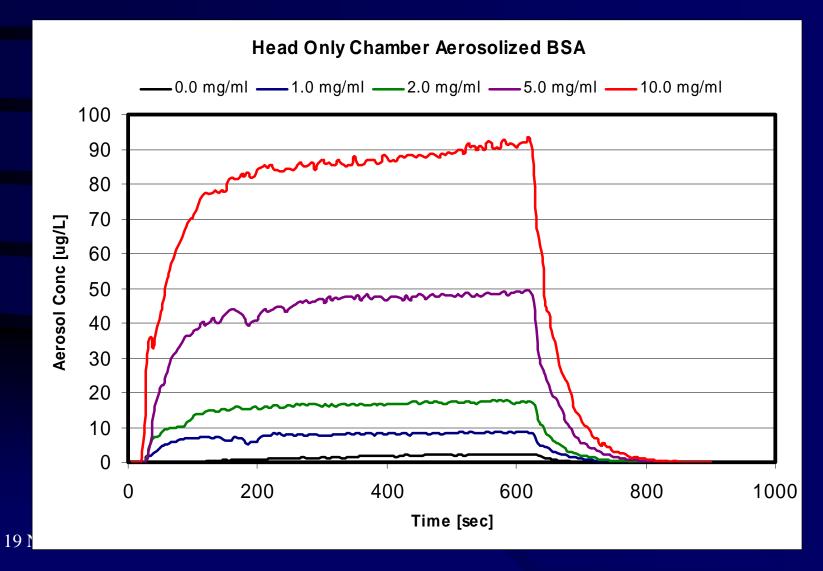


Other Biologics

- Viral and bacterial agents aerosolized in a carrier medium
- Media contain proteins
- Can media and biologic aerosol concentrations be deconvolved?

$$D(t_{\exp}) = \int_0^{t_{\exp}} R(t)C(t)dt$$

Yersinia pestis in HIB

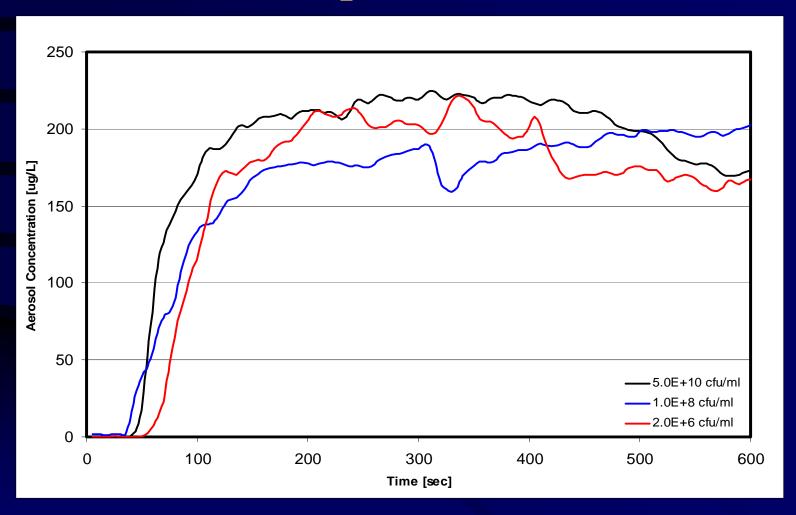

Solution: Model Concentration

- Provide C(t) mathematically
- Aerosol system then automatically compensates for R(t) variations
- Dose control through adjustment of t_{exp}

$$D(t_{\exp}) = \int_0^{t_{\exp}} R(t)C(t)dt$$

Concentration Model

Aerosol generation phase

$$C(t) = A\left(1 - e^{-kt}\right)$$

- A scaling factor related to aerosol generation efficiency
- -k time constant related to system volume and flow rates

Concentration Model

Aerosol generation phase

$$C(t) = A\left(1 - e^{-kt}\right)$$

- A scaling factor related to aerosol generation efficiency
- k time constant related to system volume and flow rates
- Chamber air wash phase

$$C(t) = C(t_{\exp})e^{-k(t-t_{\exp})}$$

 $\overline{-t_{exp}}$ – exposure time

Aerosol Model Form


- 16.3 L head-only chamber volume
- 16.0 lpm system flow rate
- Aerosol generation:

$$C(t) = A \left(1 - e^{\frac{-16.0lpm \times t}{16.3L}} \right)$$


• Chamber air wash:

$$C(t) = C(t_{\text{exp}})e^{-\frac{16.0lpm}{16.3L}(t-t_{\text{exp}})}$$

Concentration Model

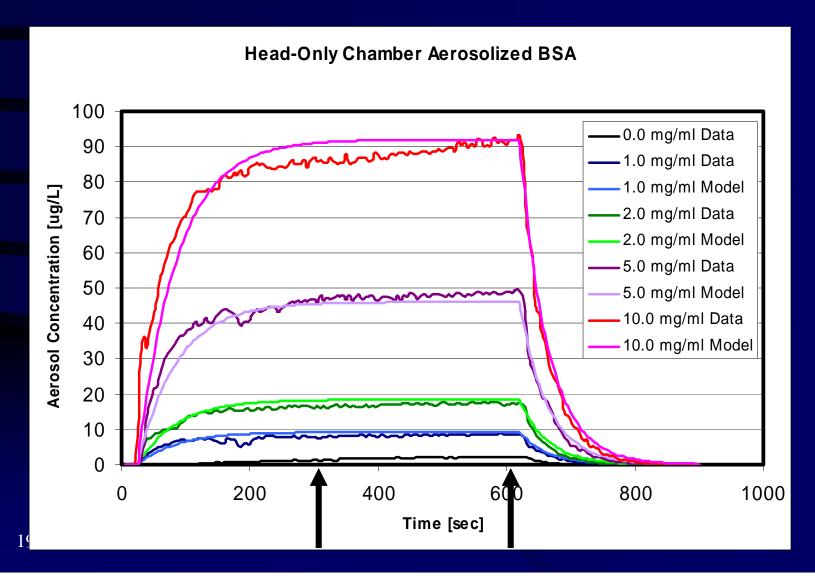

Efficiency Measurement

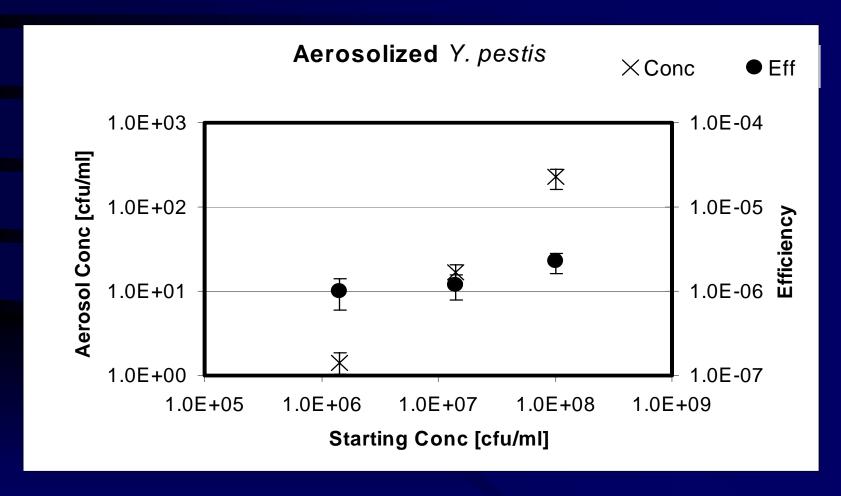

Chamber Wash Dose

- Integrated concentration over 5 minutes
- Multiplied by R(t)

$$D(t)_{wash} = \int_{t_{\exp}}^{t_{\exp}+5} C(t_{\exp}) R(t) e^{-\frac{16.0 lpm}{16.3 L}(t - t_{\exp})} dt$$

Chamber Wash Dose

- Integrated concentration over 5 minutes
- Multiplied by R(t)

$$D(t)_{wash} = \int_{t_{exp}}^{t_{exp}+5} C(t_{exp}) R(t) e^{-\frac{16.0 lpm}{16.3 L}(t - t_{exp})} dt$$

• Take R(t) as a constant, R

$$D(t)_{wash} = \frac{16.3L}{16.0lpm} RC(t_{exp})$$

$$D_{tot} = R * Eff * SC * \left\{ \left(\int_{0}^{t} 1 - e^{-\frac{16.0lpm}{16.3L}t} dt \right) + \frac{16.3L}{16.0lpm} \left(1 - e^{-\frac{16.0lpm}{16.3L}t} \right) \right\}$$

$$D_{tot} = R * Eff * SC * \left\{ \left(\int_{0}^{t} 1 - e^{-\frac{16.0lpm}{16.3L}t} dt \right) + \frac{16.3L}{16.0lpm} \left(1 - e^{-\frac{16.0lpm}{16.3L}t} \right) \right\}$$

Steady State
Concentration

$$D_{tot} = R * Eff * SC * \left\{ \left(\int_{0}^{t} 1 - e^{-\frac{16.0lpm}{16.3L}t} dt \right) + \frac{16.3L}{16.0lpm} \left(1 - e^{-\frac{16.0lpm}{16.3L}t} \right) \right\}$$

Delivered Dose

$$D_{tot} = R * Eff * SC * \left\{ \left(\int_{0}^{t} 1 - e^{-\frac{16.0lpm}{16.3L}t} dt \right) + \frac{16.3L}{16.0lpm} \left(1 - e^{-\frac{16.0lpm}{16.3L}t} \right) \right\}$$

Current Concentration


$$D_{tot} = R * Eff * SC * \left\{ \left(\int_{0}^{t} 1 - e^{-\frac{16.0lpm}{16.3L}t} dt \right) + \frac{16.3L}{16.0lpm} \left(1 - e^{-\frac{16.0lpm}{16.3L}t} \right) \right\}$$

Dose *to be*delivered during
air chamber wash


Y. pestis Aerosol Challenge

- 10 NHP
- Target dose of 75 LD_{50} (25,700 cfu)
- Minute volumes determined by plethysmography before exposure
- Starting concentration: 2e+6 cfu/ml
- Aerosol efficiency: 1.2e-6

Y. pestis Challenge

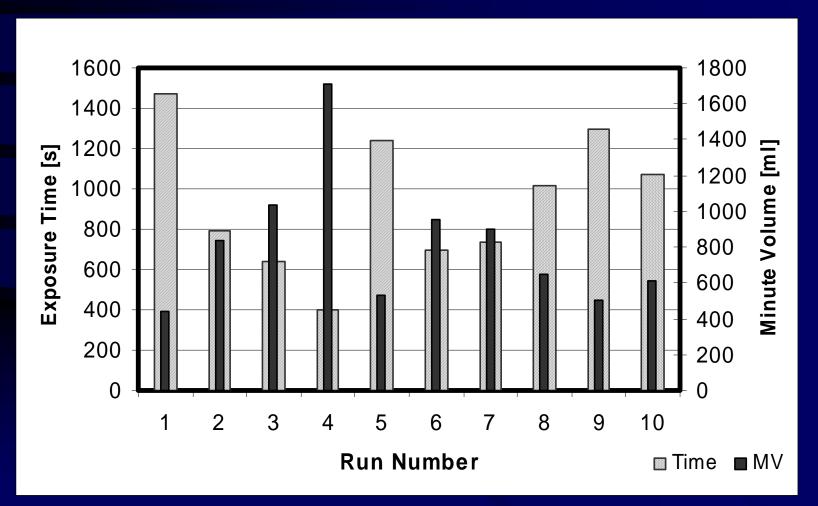

Y. pestis Challenge

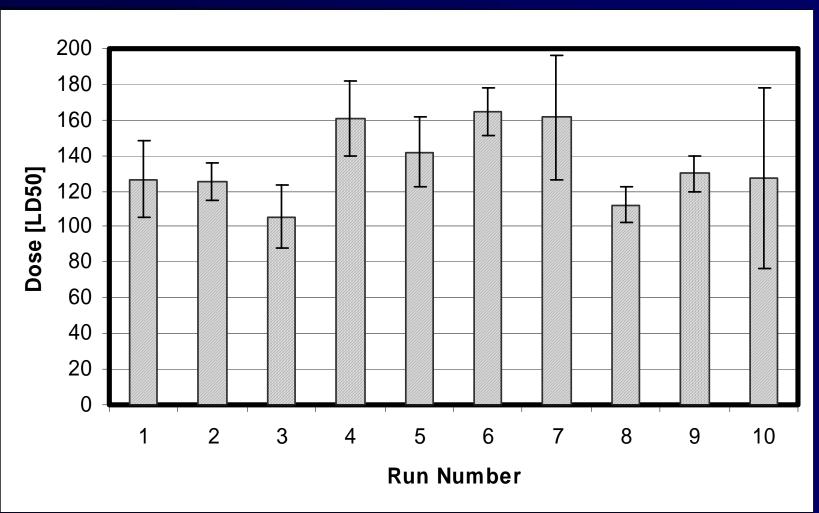

Y. pestis Challenge Results

- Using concentration model:
 - Dose range: 106 165 LD₅₀
 - $-140 \pm 20 \text{ LD}_{50}$
 - Coef of Var: 0.14
- t_{exp} range: 397 1475 sec
- Aerosol generation efficiency 2.2e-6
 - Accounts for deviation of mean dose

Comparison to Standard Exposure

- Statistics are governed by MV measurements:
 - MV range: 442 1790 ml
 - MV average: 800 ± 400 ml
 - Coef of Var: 0.5
- Dose estimation based on matching means:
 - Dose range: $74 287 \text{ LD}_{50}$
 - Dose average: 140 ± 70 LD₅₀
 - Coef of Var: 0.5

Comparison to Standard Exposure

- Statistics are governed by MV measurements:
 - MV range: 442 1790 ml
 - MV average: 800 ± 400 ml
 - Coef of Var: 0.5
- Dose estimation based on matching means:
 - Dose range: $74 287 \text{ LD}_{50}$
 - Dose average: 140 ± 70 LD₅₀
 - Coef of Var: 0.5

Three-fold reduction in C of V

Conclusions

- Improved precision of delivery of dose by aerosol to subjects with widely varying respiratory minute volumes
- Accuracy still dictated by efficiency estimates
- Three-fold reduction in coefficient of variation of dose delivery in aerosol challenge
- Isolation of biological response to aerosol challenge dose rather than to artifacts of exposure methodology

Acknowledgements

- Dr. Chad Roy
- Dr. Louise Pitt
- David Dyer
- Brittany Goodenow

Research was conducted in compliance with the Animal Welfare Act and other Federal statutes and regulations relating to animals and experiments involving animals and adheres to principles stated in the *Guide for the Care and Use of Laboratory Animals*, National Research Council, 1996. The facility where this research was conducted is fully accredited by the Association for Assessment and Accreditation of Laboratory Animal Care International.